



STARMed CaStar

Hood for CPAP therapy



Critical Care ▪ Patient Interfaces ▪ **STARMed**



STARMed CaStar CPAP hood

StarMed CaStar

The CaStar is a comfortable, versatile and lightweight CPAP patient interface designed to offer an alternative to traditional methods of CPAP delivery. The CaStar is well tolerated for extended periods and can reduce the risks associated with both endotracheal intubation and face mask use. It can be used in semi-recumbent, supine and prone positions providing an ideal option for CPAP therapy.


Technical Features

- Manufactured from phthalate-free materials
- 22M/22F connectors compliant to ISO standard
- Fitted with integrated pressure gauge to check the presence of pressure within the hood
- Comfortable and atraumatic underarm fastening straps

- Other fastening accessories are available if required
- Two versatile sealed access ports suitable for single or multiple lumen probes or catheters
- New design patient access port. Fitted with bi-directional anti-asphyxiation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated

Bibliography:

- Squadrone V, et al. Continuous Positive Airway Pressure for Treatment of Postoperative Hypoxemia. A Randomized Controlled Trial, *Jama*. 2005 Feb 2;293(5):589-95
- Squadrone V, et al. Early CPAP prevents evolution of acute lung injury in patients with hematologic malignancy, *Intensive Care Med*. 2010 Oct;36(10):1666-74
- Brambilla AM, et al. CPAP vs. oxygen therapy in severe hypoxemic respiratory failure due to pneumonia, *Intensive Care Med*. 2014 Jul;40(7):942-9
- Cosentini R, et al. Helmet continuous positive airway pressure vs oxygen therapy to improve oxygenation in community-acquired pneumonia: a randomized, controlled trial, *Chest*. 2010 Jul;138(1):114-20

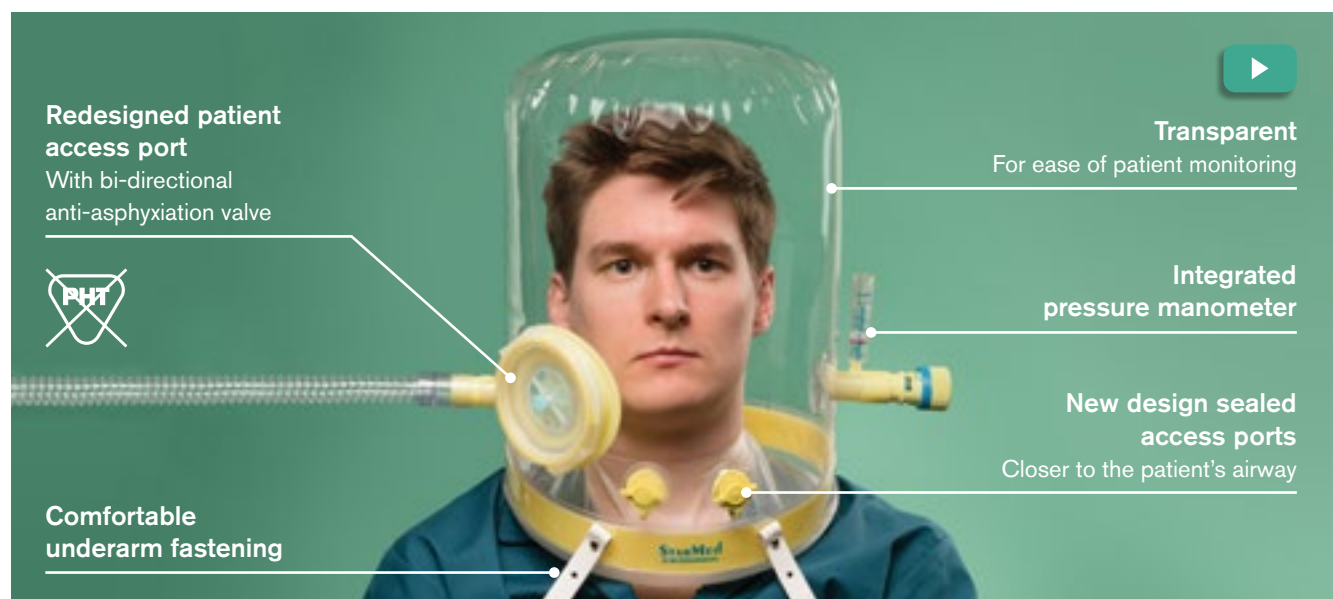

CaStar CPAP hood with bi-directional anti-asphyxiation valve and patient access port*

Code	Size	Neck circumference	Box Qty.
CP210XXL/2	Extra Extra Large	≥50cm	5
CP210XL/2	Extra Large	45-52cm	5
CP210L/2	Large	40-47cm	5
CP210M/2	Medium	34-41cm	5
CP210S/2	Small	27-34cm	5
CP210XS/2	Extra Small	17-27cm	5



www.intersurgical.nl/info/starmed

Video available at www.intersurgical.nl

Make an enquiry

*PEEP valve and tubing not included and have to be ordered separately. Please refer to the website for the range of accessories.



Intersurgical Benelux B.V. Vluchttoord 28, 5406 XP, Uden, Nederland

T: +31 (0)413 243 860 F: +31 (0)413 243 869 info@intersurgical.nl www.intersurgical.nl



Manufacturer:
Intersurgical S.p.A.
Mirandola (MO) Italy

Please think before you print

Save energy and paper. If you must print this information sheet please print it double sided.

IS1.17_NL • Issue 8 06.20