

Nasal cannula

A range of cannula for adults, pediatrics and neonates



Oxygen and Aerosol Therapy ▪ Oxygen Therapy



Nasal cannula

Intersurgical is pleased to offer a complete range of nasal cannula. This comprehensive range is entirely latex-free and includes designs suitable for use with adults, pediatrics and neonates.

Choice of options


Product choices are available to suit neonatal, pediatric and adult patients. The following prong shapes are also available:

- Straight
- Curved
- Flared
- Curved and flared


All oxygen tubing in the range is manufactured in a 'multi-channel' design, which permits oxygen flow even if the tubing is 'kinked'.

Now available with ear guards


Ear guards are an additional comfort benefit that come pre-attached to the cannula and designed for patient satisfaction.

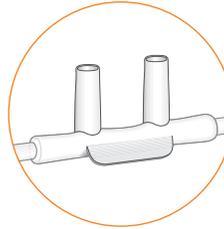


Straight Prongs

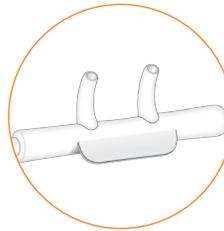

Curved prongs:
Improved anatomical shape and fit

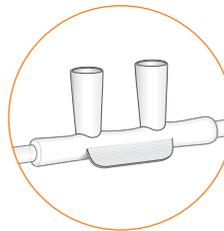

Flared prongs:
Slow down the flow of oxygen

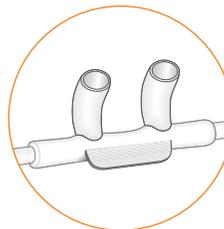

Curved and Flared prongs:
Combines the benefits of both



Curved prongs with ear guards:
Improved anatomical fit and patient comfort



Over-ear design

Secure positioning (Ear guards available for added comfort)

Non-slip toggle

Enables user to adjust for best possible fit

Multi-Channel oxygen tubing

Allows oxygen to continue to flow if the tubing becomes kinked

Clear material

Unobtrusive

Latex-free

Suitable for patients with latex allergy

Soft prongs

Minimize irritation

Location saddle

Assists in the optimum positioning of the prongs






Nasal cannula



Straight prongs

Code	Description	Length	Box Qty.
1161030	Adult, nasal cannula with straight prongs and tube	2.1m	50
1162030	Adult, nasal cannula with straight prongs and tube	7.6m	25

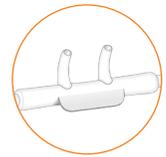


Curved prongs

Code	Description	Length	Box Qty.
1165030	Adult, nasal cannula with curved prongs and tube	2.1m	50
1165002 <i>NEW</i>	Adult, nasal cannula with curved prongs, ear guards and tube	2.1m	50
1163030	Pediatric, nasal cannula with curved prongs and tube	2.1m	50
1164030	Neonatal, nasal cannula with curved prongs and tube	2.1m	50



Flared prongs

Code	Description	Length	Box Qty.
1166030	Adult, nasal cannula with flared prongs and tube	2.1m	50

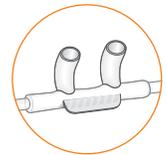

Curved and flared prongs

Code	Description	Length	Box Qty.
1167030	Adult, nasal cannula with curved/flared prongs and tube	2.1m	50