



22mm Flextube™ breathing system with integral monitoring line



Anaesthesia ▪ Anaesthetic Breathing Systems

22mm Flextube™ breathing system with integral monitoring line

This range of 22mm breathing systems are supplied with a gas monitoring line conveniently and safely positioned inside the length of the breathing system, providing a safer option.

By anchoring the monitoring line securely inside the breathing system optimal positioning is achieved for accurate and reliable gas monitoring and the need for additional luer ports in either the Y-Piece or filter/HMEF device is eliminated. This therefore avoids any possibility of accidental disconnection in use or misconnections.

The monitoring line is kept clear of the working environment, providing a tidy, less cluttered anaesthetic space and avoids the risk of occlusion or disconnection.

Combining the two products provides the convenience of a breathing system set up and ready to use and benefits the clinician by:

- providing consistent and accurate monitoring
- reducing the number of patient connections and the associated risks
- declutters the anaesthetic space around the patient

The Flextube breathing system with integral monitoring line is available in two lengths and two different configurations to suit clinical requirements.

▶ Video available at www.intersurgical.com


Convenience

One product, complete and ready to use

Accuracy

Ensures consistent and accurate CO₂ monitoring

Optimal positioning of the monitoring line at the Y-piece

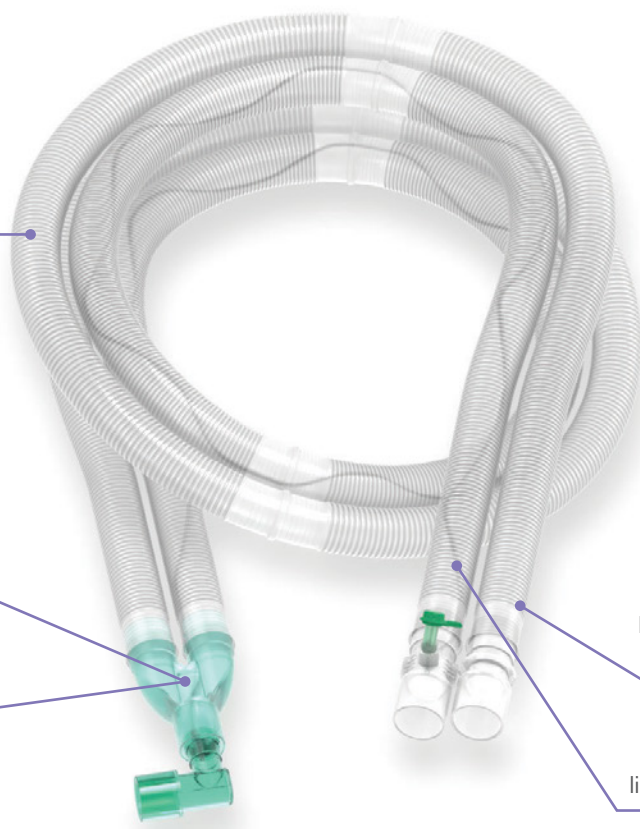

Safety

Helps to reduce the number of patient connections

Minimises risk of disconnection and occlusion



Safe and Secure

Flextube breathing systems have a unique clip feature to create an airtight fit when used with Intersurgical connectors, ensuring no leaks or disconnections

Reduces possibility of monitoring line becoming caught or entangled

Code	Description	Box Qty.
2010257	22mm Flextube breathing system with integral monitoring line, ≥1.6m	12
2010258	22mm Flextube breathing system with integral monitoring line, ≥2.4m	12
2010259	22mm Flextube breathing system with 2L bag, integral monitoring line and limb, ≥1.6m	12
2010260	22mm Flextube breathing system with 2L bag, integral monitoring line and limb, ≥2.4m	12

Make an enquiry



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
 T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016 and ISO 14001:2015

Please think before you print
 Save energy and paper.
 If you must print this information sheet please print it double sided.