



LiteScope™ laryngoscope

One adult size ▪ Fully disposable ▪ Ready for use



Airway Management ▪ Airway Accessories



Introducing the LiteScope™ laryngoscope

One adult size ▪ Fully disposable ▪ Ready for use

LiteScope is the new, one piece, single-use, fully disposable laryngoscope from Intersurgical, offering a quick and easy-to-use solution for intubation whilst maximising patient safety.

Unlike a standard laryngoscope, the innovative one-piece design of LiteScope eliminates the need to connect a separate blade and handle, thus saving valuable time in an emergency. Modelled on a Macintosh size 4 blade, the insertion technique is the same as for any standard laryngoscope with a Macintosh blade.


LiteScope is ideal for use in:

- Pre-hospital
- Emergency medicine
- Resuscitation
- Anaesthesia
- Difficult airway
- Intensive care


Infection Control

Recently published infection control guidelines have confirmed that in relation to laryngoscopes, 'current practices for decontamination and disinfection between patients are frequently ineffective, leaving residual contamination that has been implicated as a source of cross infection¹. Blades are also regularly contaminated with blood, indicating penetration of mucous membranes, and puts these items into a high-risk category.' They also confirm that, 'Ideally, a fully disposable laryngoscope, whereby the handle is also discarded with the blade, is recommended.'

As a single use, fully disposable laryngoscope, LiteScope eliminates the risks associated with the reprocessing of reusable devices, and therefore provides an ideal solution where there is a concern regarding infection control. It also avoids any reprocessing costs.

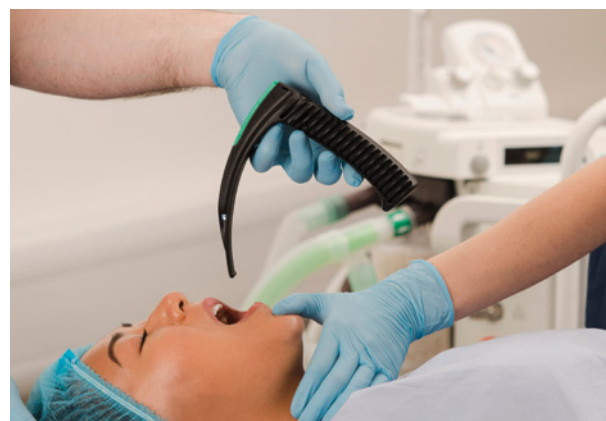

Anaesthesia

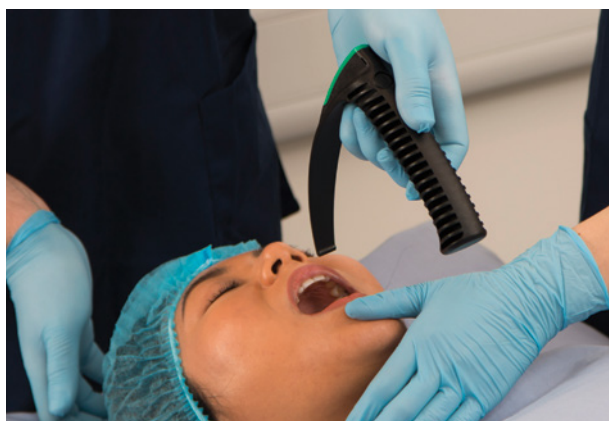

Intensive care



Emergency medicine



Features & benefits

One adult size

no need to stock multiple sizes

Fully disposable

no need to remove batteries
prior to disposal

Ready for use

No need to connect separate blades

Corrugated handle

For superior grip



On/Off switch

LiteScope is ready to use
seconds after removing
from the packaging

TPE guard

May help reduce the
risk of dental trauma
during intubation

Macintosh blade


Familiar insertion technique

Ergonomic design

Ensures LiteScope™ is
comfortable and easy to use

LED light

For optimal visualisation

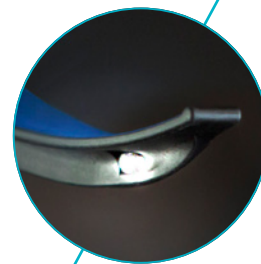

▶ Video available at www.intersurgical.com





LiteScope™ is supplied in boxes of 10 units. Each device is individually packed to ensure protection of the device, both in transit and in storage.



www.intersurgical.com/info/litescope

| Code | Description | Box Qty. |
|---------|------------------------|----------|
| 7124000 | LiteScope laryngoscope | 10 |

[Make an enquiry](#)

Intersurgical range of airway devices and accessories


The LiteScope is the most recent addition to our range of airway devices and accessories, which includes the innovative i-gel® 2nd generation supraglottic airway, the i-view™ video laryngoscope and a whole range of laryngeal masks, laryngoscope blades, tracheal tubes, oropharyngeal airways and intubating bougies and stylets. It also includes an extensive range of breathing filters, heat and moisture exchangers (HMEs) and heat and moisture exchanging filters (HMEFs). For further details regarding any of our airway products, please do not hesitate to contact us.



i-gel® supraglottic airway information sheet




i-view™ video laryngoscope information sheet

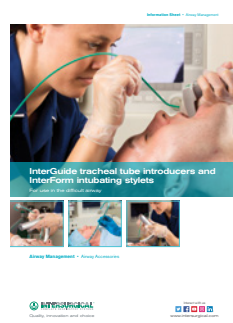

InterGuide tracheal tube introducers and InterForm intubating stylets information sheet



Safe and effective breathing filters information sheet

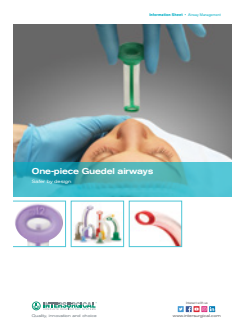

One-piece Guedel airways information sheet



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016 and ISO 14001:2015

Please think before you print
Save energy and paper.
If you must print this information sheet please print it double sided.

IS6.17 • Issue 1 03.21