



ClearFlex™ Plus anaesthetic face mask





ClearFlex™ Plus anaesthetic face mask

The ClearFlex Plus anaesthetic face masks are manufactured from silicone with an integrated plastic shell. This design helps to facilitate a comfortable grip and provide mask stability. Due to the anatomical shape, the ClearFlex Plus mask can provide a leak-free seal with minimum pressure and the ergonomically design allows an optimised grip via C-grip. The ClearFlex Plus mask is available in 5 sizes and can be reprocessed up to 100 times (autoclavable at 134°C).

Safe and secure connection

The 22 mm tapered connection allows connection to either resuscitator bag or breathing system

Reinforced design

The over moulded reinforced body helps to provide mask stability and supports a comfortable grip

UDI code and serial number

Allowing product traceability and tracking of product usage



Clear mask body


Provides clear visibility for confirmation that the patient is well perfused

Clear sizing and colour coding

Colour coded size numbers facilitates easy size selection

Soft silicone cuff

Seals to a variety of face shapes, providing a comfortable and effective seal

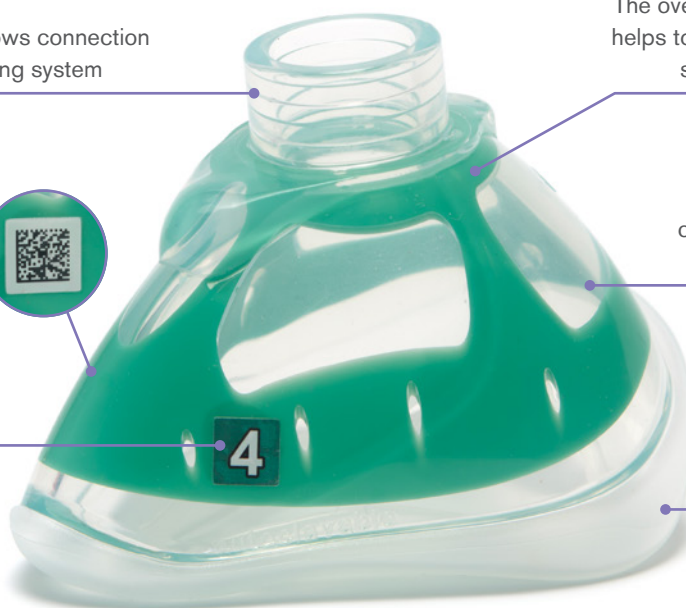

Code	Description	Colour	Box Qty.
8749666	ClearFlex plus silicone anaesthetic face mask, size 6, extra large adult, 22F		1
8749555	ClearFlex plus silicone anaesthetic face mask, size 5, large adult, 22F		1
8749444	ClearFlex plus silicone anaesthetic face mask, size 4, medium adult, 22F		1
8749333	ClearFlex plus silicone anaesthetic face mask, size 3, small adult, 22F		1
8749222	ClearFlex plus silicone anaesthetic face mask, size 2, paediatric, 22F		1

Note: Not all products are available in all markets, please contact your local office for details.

IS5.47 • Issue 1 06.26



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to
ISO 14001:2015, ISO 9001:2015, ISO 13485:2016 and MDSAP

UK • Australia • Benelux • Brazil • Canada • China • Colombia • Czech Republic • Denmark • France • Germany • India • Ireland • Italy
Japan • Lithuania • Malaysia • Mexico • Norway • Philippines • Portugal • Russia • South Africa • Spain • Sweden • Taiwan • Turkey • USA