




# STARMed Castar Up

Openable CPAP hood with inflatable neck cushion

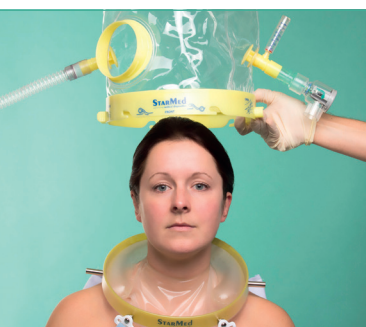

**Critical Care** ▪ Non-invasive ventilation



# STARMed CaStar Up – openable CPAP hood

## StarMed CaStar Up

Openable Castar Up is a device recommended for the treatment of a number of respiratory conditions, including:


- Cardiogenic pulmonary oedema
- Non cardiogenic pulmonary oedema (ALI-ARDS)
- Thoracic trauma
- Atelectasis
- Post operative hypoxemia
- Asthma
- COPD exacerbation
- Pre and post extubation respiratory support

## Advantages

- Allows rapid access to the patient's face during nursing procedures or in case of emergency
- One operator can easily remove the upper part of the hood and have full access to the patient
- Probes and catheters can remain in place, on the patient, connected to the lower part of the hood
- Well tolerated for extended periods
- Prevents the risk of pressure sores from face masks

## Technical Features

- 22M/22F connectors compliant to ISO standard
- Fitted with integrated pressure gauge to check the presence of pressure within the hood
- Internal volume 15 litres – weight 690g
- Inflatable cushion for added patient comfort
- Comfortable underarm fastening straps
- Other fastening accessories are available if required
- Two new larger access ports on the collar for probes or catheters of 3.5-7.0mm Ø
- Patient access port with screw cap. Fitted with bidirectional anti-suffocation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated
- Additional models available with a built-in inflatable neck cushion or without patient access port

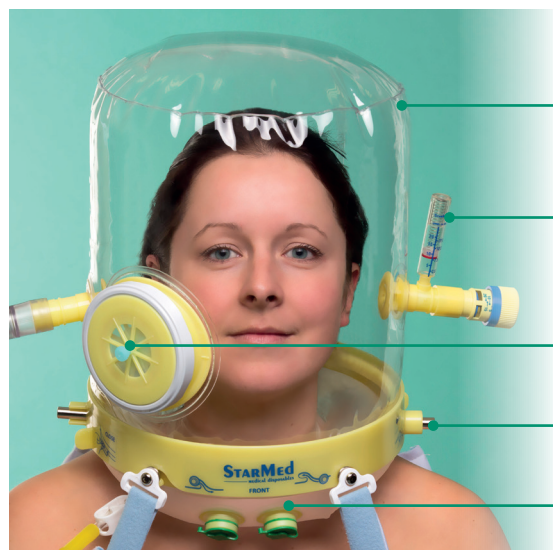

## Features and benefits

### Transparent

For ease of patient monitoring

### Integrated pressure manometer

### Patient access port with anti-suffocation valve

### Easy removal for rapid patient access

### New large dimension accesses on the collar

Ideal for single and multiple lumen probes or catheters

## CaStar Up CPAP hood with anti-suffocation valve and patient access port

Code	Description	Box Qty.
CA011XXL/2	Extra Extra Large ≥50cm	5
CA011XL/2	Extra Large 45-52cm	5
CA011L/2	Large 40-47cm	5
CA011M/2	Medium 34-41cm	5
CA011S/2	Small 27-34cm	5
CA011XS/2	Extra Small 17-27cm	5

Single packed. Minimum purchasing package: 5 units. Custom-packs are supplied on specific request of the customer.



CE certified product

IS1.31\_ROW • Issue 2 02.14



Intersurgical Ltd,  
Crane House, Molly Millars Lane,  
Wokingham, Berkshire, RG41 2RZ, UK  
T: +44 (0)118 9656 300 F: +44 (0)118 9656 356  
info@intersurgical.com www.intersurgical.com

Deutschland  
info@intersurgical.de

France  
info@intersurgical.fr

Italia  
info@intersurgical.it

España  
info@intersurgical-es.com

Portugal  
info@intersurgical.pt

Nederland  
info@intersurgical.nl

Lietuva  
info@intersurgical.lt

Россия  
info@intersurgical.ru

Česká Republika  
info@intersurgical.cz

Taiwan  
info@intersurgical.com.tw

South Africa  
info@intersurgical.co.za

Philippines  
info@intersurgical.ph

Japan  
info@intersurgical.co.jp

USA  
support@intersurgicalinc.com

China  
info@intersurgical-cn.com

