




Transfer tube safety valve

For Anaesthesia Gas Scavenging System (AGSS)

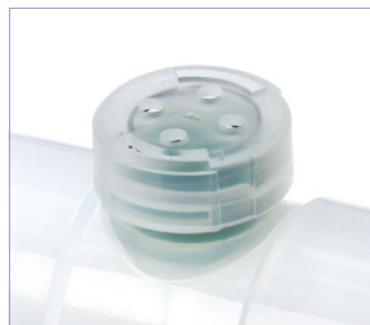

Anaesthesia ▪ Breathing System Accessories

Scavenging safety valves

Due to the potential danger of excessive pressure caused by a blockage in the AGSS transfer tube, the MHRA (Medicines and Healthcare products Regulatory Agency UK) issued a Medical Device Alert (MDA/2010/021) calling for the inclusion of a safety pressure relief valve.

Subsequently the European and International standard EN ISO 80601-2-13 Anaesthetic workstation in section 201.103.3.2.1 has been updated to reflect this alert.

Wherever there is a 30mm port on the anaesthetic breathing system, for example on the exhaust outlet of the APL valve, or any gas exhaust port of the anaesthetic machine, ventilator or absorber, one of these valves must be fitted as a safeguard in case of an occlusion in the transfer tube.



Video available at www.intersurgical.co.uk



Code	Description	Box Qty.
1972000	Anaesthesia gas scavenging pressure relief valve 30F/30F with 4.2M 30mm Flextube™	5
1972001	Anaesthesia gas scavenging pressure relief valve 30F/22M with 4.2M 22mm Flextube	5

Make an enquiry

Scavenging connectors

A range of connectors to use with our 22mm and 30mm scavenging tube, allowing the tube to be connected to the breathing system exhaust valve, or anaesthesia machine exhaust valve and the AGSS device.


The ISO standard is 30mm, and the ANSI standard is 19mm. Historically 22mm has also been used.

				
Code 1970000	Code 1970001	Code 1971000	Code 1971001	Code 1975000
Box Qty. 35	Box Qty. 75	Box Qty. 25	Box Qty. 75	Box Qty. 50

Make an enquiry

Scavenging tube


Tinted blue to differentiate it from the anaesthesia breathing system tubing used in operating theatres, the 30mm scavenging tube conforms to ISO standards for gas scavenging.

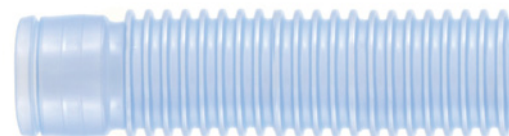

Video available at www.intersurgical.com



 illustrated

Code	Description	Box Qty.
 1574066	22mm Flextube blue scavenging tube – 50m roll cuffed at 0.4m	1
1573000	22mm Flextube – 50m roll cuffed at 0.15m	1
1574000	22mm Flextube – 50m roll cuffed at 0.4m	1
1570000	30mm Flextube blue scavenging tube – 30m roll cuffed at 1.6m	1

Make an enquiry



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK

T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016 and ISO 14001:2015

Please think before you print
Save energy and paper.
If you must print this information sheet please print it double sided.

IS 5.27 • Issue 4 07.19