

Explorer masks

Self-sealing valves

To facilitate easy access



Explorer masks

Explorer masks for single patient use to prevent the possibility of cross-infection.

The Intersurgical Explorer masks are designed to allow patient ventilation and access simultaneously during procedures.

There are three one-way valves positioned on the mask (with the exception of size 1 where there is only one valve), which provide easy access from either side of the patient.

The masks are primarily used during procedures such as fiber optic intubation, endoscopy, bronchoscopy and T.E.E. procedures. Ventilation is delivered through the central 22mm connection port.

Features and benefits

Color-coded hook ring

For ease of identification and integration with other Intersurgical airway products

Anatomical shape

For a precise fit

4 sizes available

From Adult to Infant






Soft cushion

For patient comfort and air-tight fit

Self-sealing valves

To facilitate easy access

Ordering information

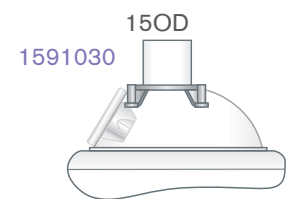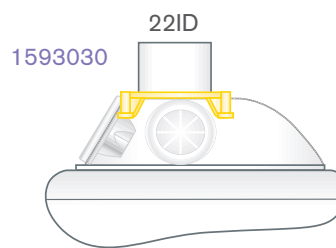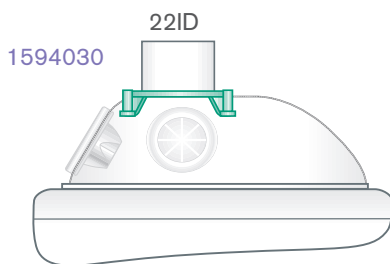

| | | | | |
|---------------------------|--------|---|----------|-------------|
| Code: 1591030 | | | | Box qty.50 |
| Infant Explorer Mask | size 1 | ■ | 1 valve | |
| Code: 1592030 | | | | Box qty. 50 |
| Pediatric Explorer Mask | size 2 | □ | 3 valves | |
| Code: 1593030 | | | | Box qty. 50 |
| Small Adult Explorer Mask | size 3 | ■ | 3 valves | |
| Code: 1594030 | | | | Box qty. 50 |
| Adult Explorer Mask | size 4 | ■ | 3 valves | |