

Nasal Cannula

A range of cannula for adults, pediatrics and neonates



Nasal Cannula

Intersurgical is pleased to offer a complete range of nasal cannula. This comprehensive range is entirely latex-free and includes designs suitable for use with adults, pediatrics and neonates.

Choice of options


Product choices are available to suit neonatal, pediatric and adult patients. The following prong shapes are also available:

- Curved
- Flared
- Curved/Flared
- Straight


All oxygen tubing in the range is manufactured in a 'multi-channel' design, which permits oxygen flow even if the tubing is 'kinked'.

Now available with ear guards


Ear guards, an additional comfort benefit, are pre-attached to the cannula and designed for patient satisfaction.

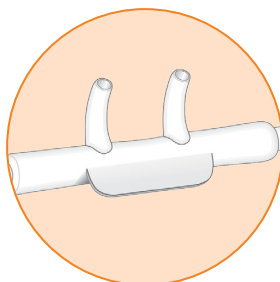


Curved prongs:
Improved anatomical fit

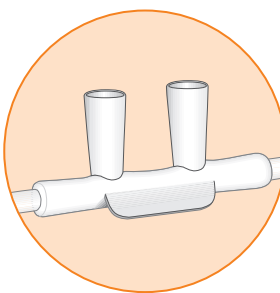

Flared prongs:
Slows down the flow of oxygen

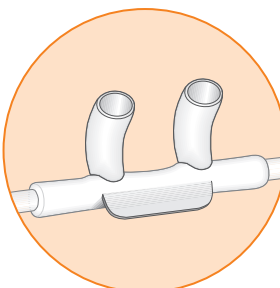

Curved / Flared prongs:
Combines the benefits of both the above designs

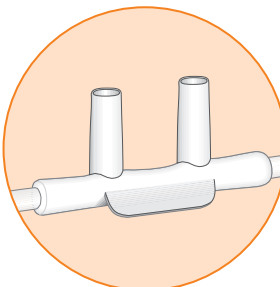

Straight prongs

Features and Benefits

Multi-Channel oxygen tubing

Allows oxygen to continue to flow if the tubing becomes kinked

Latex free

Suitable for patients with latex allergy

Over-ear design

Secure positioning
(Ear guards available for added comfort)

Non-slip toggle

Enables user to adjust for best possible fit

Clear material

Unobtrusive

Soft prongs

Minimize irritation

Location saddle




Assists in the optimum positioning of the prongs








Adult headset version

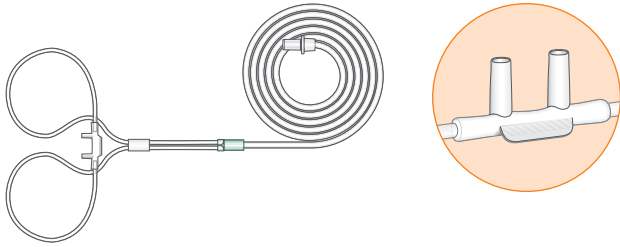

Ideal for domiciliary use

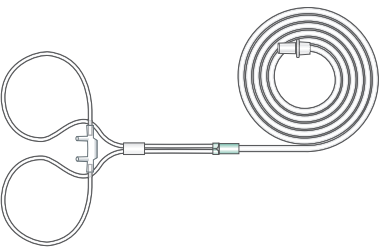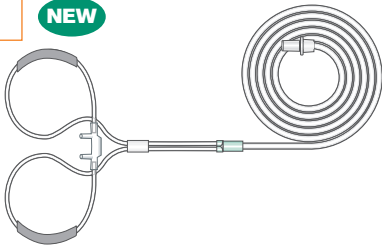


Ordering information

<div>1161030</div> 		
1161030	Adult Straight Prong with tube, 7' length	50 
1162030	Adult Straight Prong with tube, 25' length	25 

<div>1165030</div>  <div>1165002 NEW</div>  		
1165002	Adult Curved Prong with ear guards and tube, 2.1m length	NEW 50 
1165030	Adult Curved Prong with tube, 7' length	50 
1164030	Neonatal Curved Prong with tube, 7' length	50 
1163030	Pediatric Curved Prong with tube, 7' length	50 
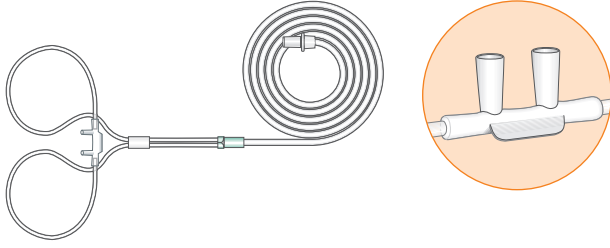
<div>1166030</div> 		
1166030	Adult Flared Prong with tube, 7' length	50 
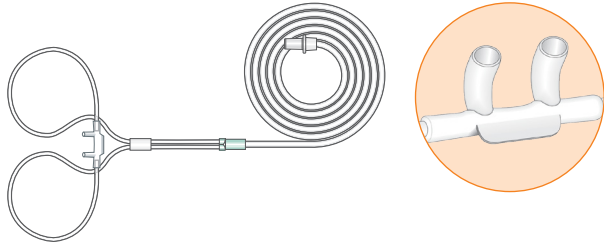
<div>1167030</div> 		
1167030	Adult Curved / Flared Prong with tube, 7' length	50 